

# BALMORAL PLUS ACTIVE MATTRESS



**Product Code:**

HAMBALM004

**Description:**

- The mattress is made up of 27 individual cells; 23 alternating air cells and 4 static pillow section cells
- The mattress has an integrated air pocket base
- Cells are made of pure PU
- Pillow section for patient stability
- All cells are interchangeable for ease of repair and cleaning
- 4-way stretch welded PU cover which is vapour permeable and water resistant and has a Crib 5 fire rating
- The Balmoral Plus has been designed with a three-part static air cell mat that can be deflated in the event of cardiac arrest by opening the CPR valve.
- Compatible with Harvest Healthcare's pumps 250 and 3

**Mattress Dimensions:**

2000 x 900 x 150 mm

**Weight:**

10kg

**Max. Patient Weight:**

Harvest 250 Pump: 31 Stone (200kg)  
Harvest 3 Pump: 40 Stone (254kg)

**Pressure Relief Risk:**

Very High

**Pump Technical Specification:**

Please see separate pump specification sheet

**Cleaning Guidance:**

Cover should be wiped down with hot water and mild detergent. If heavily soiled, the cover may be machine washed up to 60-95°C and tumble dried. For disinfection use sodium hypochlorite solution max 1 000ppm. Use no bleach, biological or phenolic-based cleaners. Dry completely before storing.

**Warranty**

In terms of Harvest Healthcare Guarantee and Warranties the Harvest Balmoral Plus mattress is guaranteed for three years from date of purchase.

**Exclusions:**

- Damage to fabric covers through incorrect use and penetration by sharp instruments is excluded.

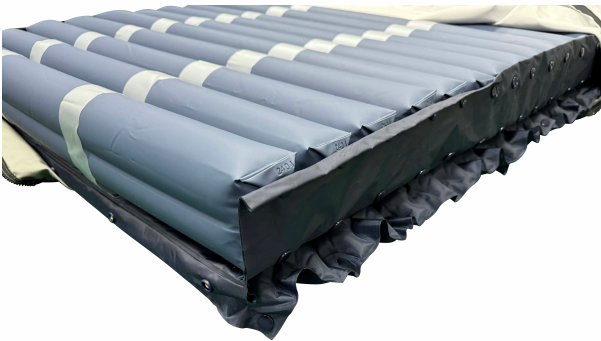
**Supplier Contact Details:**

t: +44(0)1709 377172

e: [sales@harvesthealthcare.co.uk](mailto:sales@harvesthealthcare.co.uk)

w: [www.harvesthealthcare.co.uk](http://www.harvesthealthcare.co.uk)

Sheaf house, Bradmarsh Way, Bradmarsh Business Park, Rotherham UK, S60

**Balmoral Plus cells****Important Safety Notice**

Harvest Healthcare Ltd supplies equipment and accessories to cover a range of medical conditions. Specific choice of equipment is the sole responsibility of the clinical practitioner or carer overseeing individual patient care. Harvest Healthcare Ltd supplies a range of products such as foam mattresses, active air mattresses, table protection mats and bed rail protection bumpers. These products are usually covered in vapour permeable PU fabric that could cause risk of suffocation in the event of prolonged table contact. Mattress risk categories are to be used as a guide only. RSK assessments specific to user must be carried out prior to the use of all products. Pre-UKCIs and data18 content may undergo development and may therefore deviate from the details given.

**Version 1: June 2024**